

GMI GABIPOORT IMPLANTABLE CATHETER

1500 - VASCULAR LINE



ACCURACY

- Trilobed septum that facilitates the correct location of the puncture.

VERSATILITY

- Different combinations of sizes between chambers and tubes.
- Low profile and high flow catheter.

RELIABILITY


- Developed according to ISO 10555-6.

SAFETY

- Excellent visualization in imaging exams.

SAFETY CONNECTION


- Ergonomic piece.
- Safe and easy locking that ensures effective use.

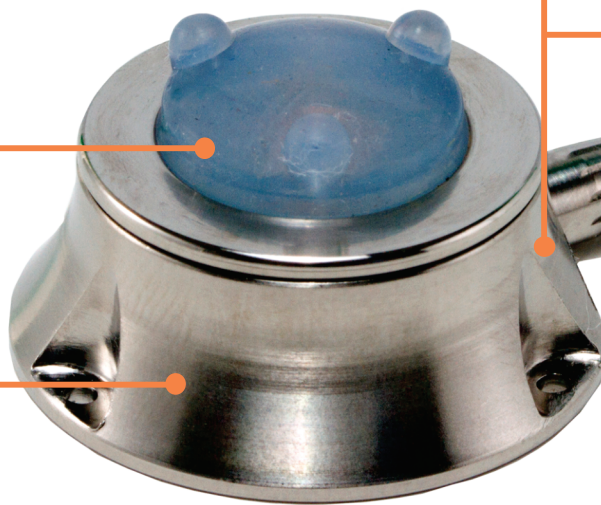

GABI
PORT
PROVIDES A
COMPLETE
AND
ADAPTABLE
SOLUTION
FOR THE
DIFFERENT
PATIENTS
PROFILES
THAT
NEED
PROLONGED
VENOUS
ACCESS.

DIVERSITY

-Different sizes for better adaptation and safety in patients.

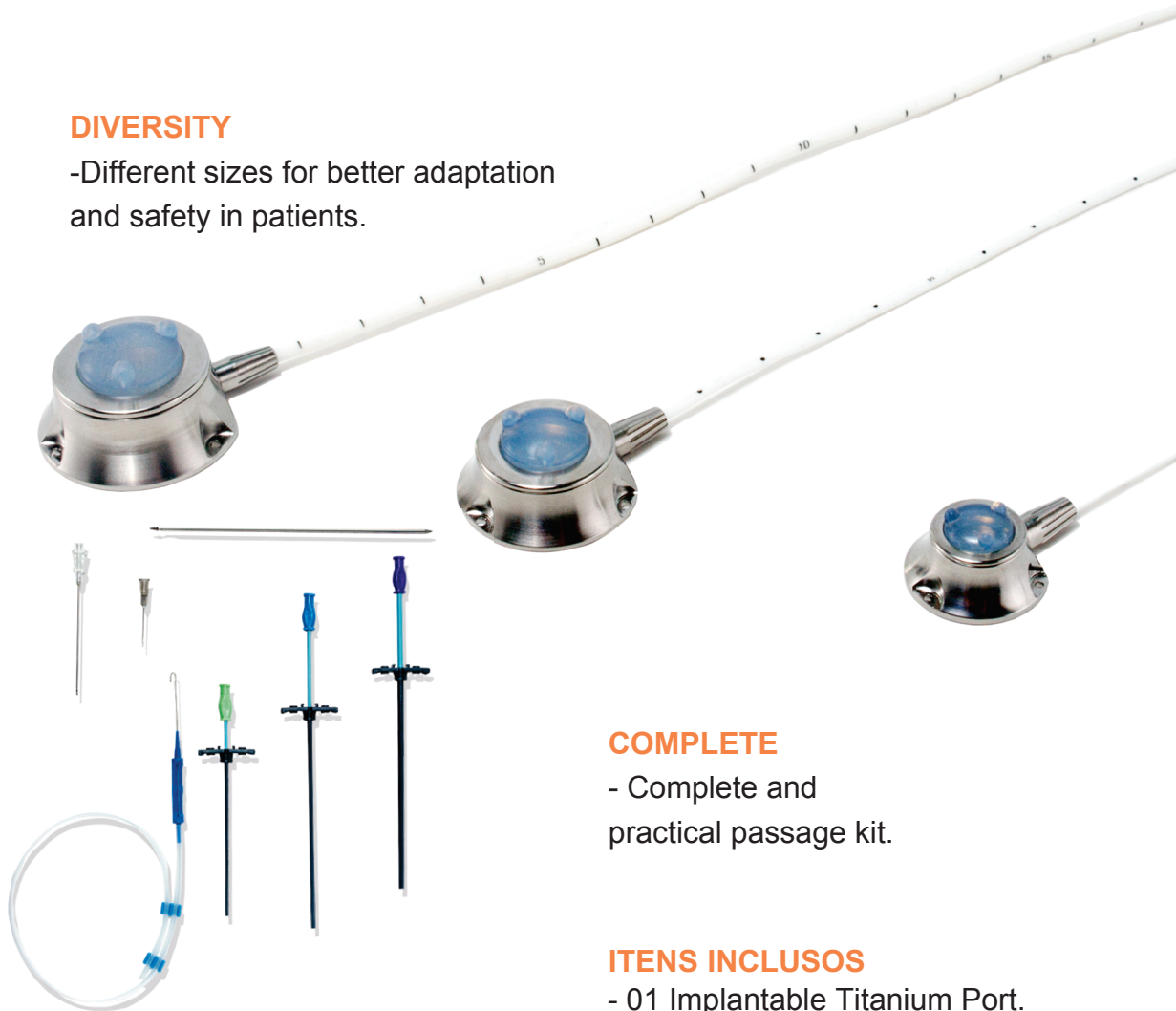

COMPLETE

- Complete and practical passage kit.

ITENS INCLUSOS

- 01 Implantable Titanium Port.
- 01 Silicone Catheter.
- 01 Straight Huber Needle.
- 01 Introducer Needle.
- Stainless Steel Guide Chuck with Door J Type and Protective Cover.
- 01 Tunneler.
- 01 Peel Away Dilator.
- 01 10ml syringe with *luer-lock*.

INDICATION

Prolonged vascular access for infusion of intravenous therapy, chemotherapy, antibiotic therapy, parenteral nutrition and blood products. It suits the need for patient therapy, due to the diversity of sizes.

	1500 GMI FULLY IMPLANTABLE CATHETER ANVISA Registration No. 80423540030					
	CHAMBER			SEPTUM		
Code	INDICATION	DIAMETER x HEIGHT (mm)	PRIME VOLUME (ml)	EXTERNAL DIAMETRO (mm)	THICKNESS (mm)	Catheter (fr x cm)
1500-10-04	Lac/Ped	19,0x8,4	0,18	9,0	4,1	4,0x60
1500-10-05	Lac/Ped	19,0x8,4	0,18	9,0	4,1	5,0x60
1500-20-05	Ped/Adult	25,0x9,3	0,40	12,0	4,5	5,0x60
1500-20-07	Ped/Adult	25,0x9,3	0,40	12,0	4,5	7,0x60
1500-20-09	Ped/Adult	25,0x9,3	0,60	12,0	4,5	9,0x60
1500-30-07	Adult	28,5x11,7	0,60	14,5	6,6	7,0x60
1500-30-09	Adult	28,5x11,7	0,60	14,5	6,6	9,0x60